

WATER KITS SUPPLY

Téléport 5 - 62590 Juillan - France

Tél: +33(0) 5 62 95 17 94 Fax: +33(0) 5 62 95 34 27 Email: contact@water-kits.fr Web: www.water-kits.fr

Coliform/E.coli Test Kit

- Early identification of bacterial contamination
- Simple one step procedure
- Coliforms and E.coli in one test
- 100ml sample (regulatory reporting)
- Detects one CFU/100ml
- Short incubation and detection time (18-24hours)
- No Sodium Thiosulphate needed
- 2 year shelf-life at room temperature

The test is available in pre-packed bottle form – where powder is loose in a 100ml bottle






Color change (Green-Blue)

= Coliform positive



A Mini Fluorescent lamp is required for the determination of E.coli.



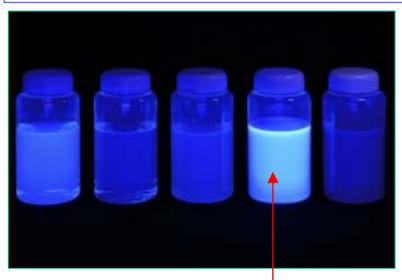

Observation of fluorescence might = *E. coli* positive

CF80	Coliform/E.coli bottles (100ml x 80 bottles)
INC-01	Compact field Incubator
BI-K27718	Mini Fluorescent lamp